Cedis Profi UVC CLEAN

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Internet: www.egger-labor.com

Cedis Profi UVC CLEAN

The professional device for hygienic certainty in your acoustician shop.

User instructions

**Read these instructions carefully prior to use
and keep them in a safe place!**

General product information:

Sterilizing surfaces with UVC light eliminates over 99% of viruses and bacteria. The UVC lighting technology destroys microorganisms, germs, bacteria and fungi. The ring-shaped UVC lamp (360°) together with the high-gloss metallized chamber makes this cleaning possible on the entire surface.


Cedis Profi UVC CLEAN was especially developed for the use as described above and is characterized by:

- simple operation
- easy access due to hinged lid
- 3 and 5 minute timer function with automatic safety switch-off upon opening lid
- blue LED indicators from left to right display the operation time i.e. the end of the sterilizing process

Prior to first operation:

- Place the device on a dry and sturdy surface before operating
- Do not submerge the device in water
- Only use the device for cleaning earmolds, BTE and ITE hearing aids
- Do not use the device on animals! Only use the device for cleaning objects!
- Keep away from heat sources
- Keep away from children

Product description:

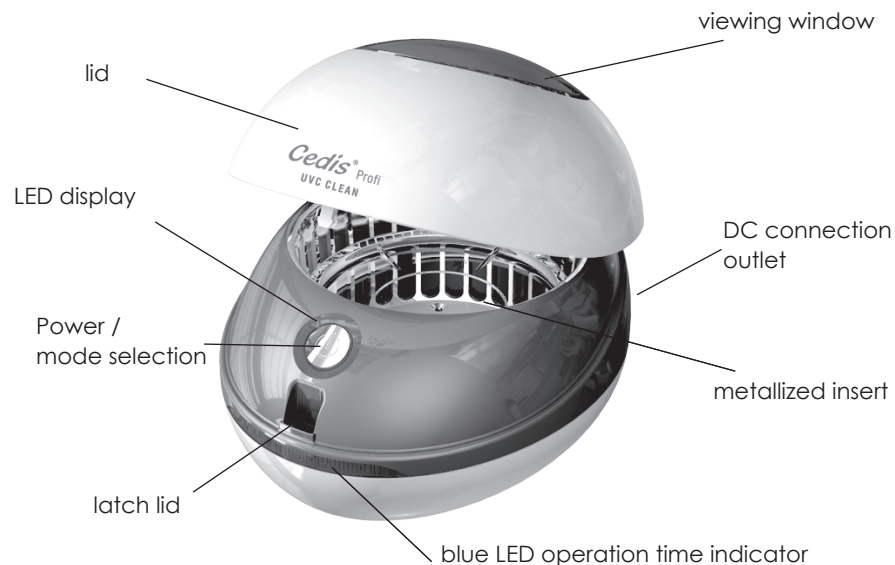Product specifications:

UVC light tube: 4 W ($\pm 15\%$) DC

UVC wavelength: 253.7 nm (ozone free)

Power unit: Input: 100~240 V – 50/60 Hz
Output: 6 V

LED display: Blue LED light to indicate 3 minutes /
Green LED light to indicate 5 minutes

Outer dimensions: W 143 x D 167 x H 108 mm

Inner dimensions: Ø 105 mm, H 55 mm

Scope of delivery:

Cedis Profi UVC CLEAN

Power supply unit

User instructions

User instructions:

1. Insert the jack plug of the power supply unit (12 V), included in delivery, in the connector at the back of the device. Plug the power supply unit into a suitable electrical outlet - your Cedis Profi UVC CLEAN is now ready for operation.
2. Gently press down on the lid and it will spring up automatically. Lift the lid.
3. Open the battery compartment of BTE and ITE hearing aids or switch off your hearing system.
4. Place the earmolds, BTE or ITE hearing aids into the metallized insert of the device.
5. Press the power button to select the operation time:
TIP: Greatest effectiveness of more than 99% is achieved by sterilizing for 5 minutes.
 - Press the button once to activate sterilizing for 3 minutes (LED display will turn BLUE).
 - Press the button twice to activate sterilizing for 5 minutes (LED display will turn GREEN).
 - Press the button three times to deactivate the device (LED display will turn off).
6. Close the lid until you hear a "click". Due to the securing of the lid, the device is activated and the sterilizing process is started.
7. The operation time is displayed by 4 blue LED indicators that gradually turn off from left to right during sterilizing. (3-minute mode 45 seconds per LED indicator, 5-minute mode 75 seconds per LED indicator).
8. After the dimming of the four blue LED indicators, the device turns off automatically. After opening the lid, the earmolds, BTE or ITE hearing aids can be removed.

The device is equipped with an automatic safety switch-off function that is activated when opening the lid during operation. Upon closing the lid, the selected operation mode will start over.

Important safety warnings:

- Do not clean your device using solvents
- Do not clean your device under running water
- Do not pull the cable to detach the power supply unit from the electrical socket
- If the cable is damaged, it must be replaced immediately
- Repairs on the device are only to be executed by specialized personnel

Warranty:

The warranty of the device is valid for 12 months.
Please keep the receipt.

Attention:

In case of improper handling of the Cedis Profi UVC CLEAN all warranty and liability claims will become extinct.

Customer service:

In case of damages and defects please contact your distributor or

egger Otoplastik + Labortechnik GmbH
Aybühlweg 59
87439 Kempten/Germany
Phone international: +49 831 58113-60
Fax international: +49 831 58113-14
E-mail: sales@egger-labor.de
Internet: www.egger-labor.com

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